

Connect with a NewLeaf Home Medical Advisor at 888.739.4330 for personal assistance or to order your very own Invacare micorAIR mattress today!



Shop at NewLeafHomeMedical.com

INVACARE® MICROAIR® MA600
THERAPEUTIC SUPPORT SURFACE



Yes, you can.®

microAIR®

Invacare® microAIR® MA600

Alternating Pressure Low Air Loss Mattress System

Medium- High Risk
Stage I-IV

3-in-1 Alternating Pressure System with
Low Air Loss aids in the Prevention and
Treatment of Pressure Injuries



Highly vapor permeable and oversized pliable quilted nylon top cover provides low shear, friction and moisture protection



Cable management system keeps cords away from the ground to help prevent trip hazards




Nylon non-slip base cover
helps reduce shear and friction during bed articulation




Features & Benefits

8 adjustable comfort settings depending on the patient's therapy requirement

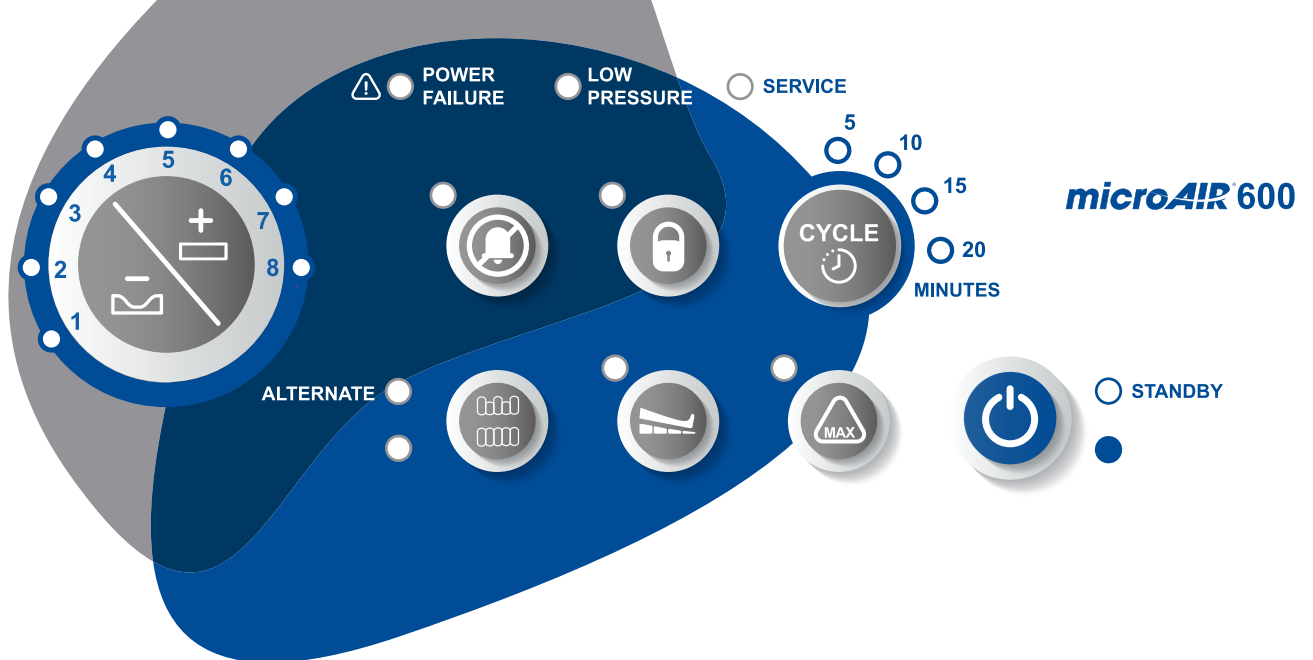

Advanced 3-in-1 alternating function inflates and deflates cells in a 3:1 cycle, providing active pressure relief by supporting 2/3 of the body at any given time

Max Inflate Mode provides quick maximum inflation within seconds to aid with transfers and nursing procedures

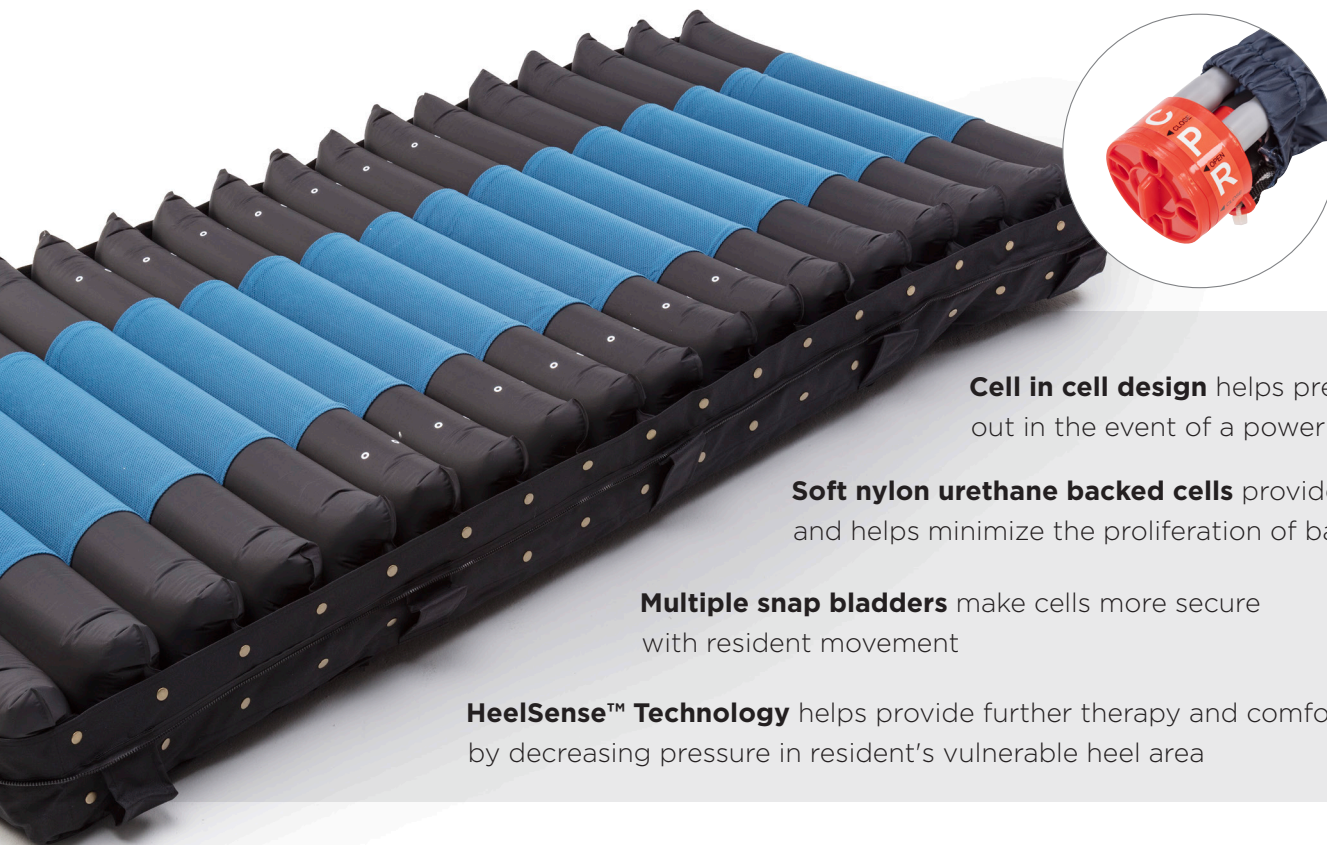

Cell in cell design helps prevent bottoming out in the event of a power failure

Soft nylon urethane backed cells provide easy cleaning and helps minimize the proliferation of bacterial growth

Multiple snap bladders make cells more secure with resident movement

HeelSense™ Technology helps provide further therapy and comfort by decreasing pressure in resident's vulnerable heel area

Request a Quote or Place an Order Today!



Contact your Inside Sales Representative

or



Contact Customer Service at 800.333.6900

Invacare Corporation

www.invacare.com

USA

One Invacare Way
Elyria, Ohio
44035-4190
(800) 333-6900

Canada

570 Matheson Blvd. E.,
Unit 8
Mississauga, Ontario
L4Z 4G4 Canada
(800) 668-5324


©2019

Invacare Corporation.
All rights reserved.
Trademarks are identified by the symbols™ and ©. All trademarks are owned by or licensed to Invacare Corporation unless otherwise noted. Specifications are subject to change without notification.
Form No. 19-024
190203



Yes, you can.®

microAIR MA600 Alternating Pressure Low Air Loss Mattress System

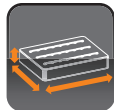

Model #	Dimensions	Weight Capacity	Flow Rate (Pump)	Therapy Type	Alarms	Limited Warranty
MA600	80" L 36" W 8" H	500 lb.	10 LPM	Alternating Pressure Low Air Loss	Audio & Visual	1 Year (Cover/Mattress) 1 Year (Power Unit)
MA600M	microAIR MA600 Mattress only					
MA600P	microAIR MA600 Pump only					