

Connect with a NewLeaf Home Medical Advisor at 888.739.4330 for personal assistance
or to order your very own Invacare micorAIR mattress today!



Shop at [NewLeafHomeMedical.com](https://www.NewLeafHomeMedical.com)

INVACARE® MICROAIR® MA800
THERAPEUTIC SUPPORT SURFACE



microAIR®

Invacare® microAIR® MA800

Alternating Pressure Low Air Loss Mattress System

**Medium - High Risk
Stage I-IV**

3-in-1 Alternating Pressure combined with Pulsation and True Low Air is Loss ideal for situations where maceration and moisture management are an issue



Highly vapor permeable and oversized pliable quilted nylon top cover provides low shear, friction and moisture protection



Cable management system keeps cords away from the ground to help prevent trip hazards



Nylon non-slip base cover helps reduce shear and friction during bed articulation





Features & Benefits

10 adjustable comfort settings depending on the patient's therapy requirement

60 second pulsation feature helps improve capillary blood flow by continually readjusting to a level surface pressure

Auto Firm Mode provides quick maximum inflation within seconds to aid with transfers and nursing procedures

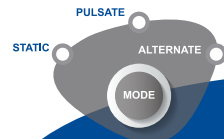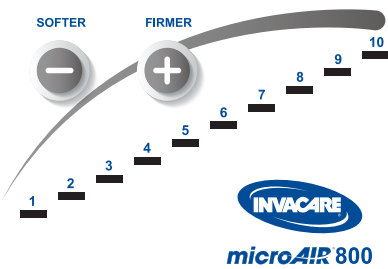

AUTO FIRM

LOCK OUT

ALARM RESET

STANDBY

OPERATE

Advanced 3-in-1 alternating function inflates and deflates cells in a 3:1 cycle, providing active pressure relief by supporting 2/3 of the body at any given time

Customizable cycle times to meet patient's varying needs

True Low Air Loss Therapy provides extra comfort for moisture and temperature control



Foam base helps prevent bottoming out in the event of a power failure

Soft nylon urethane backed cells provide easy cleaning and helps minimize the proliferation of bacterial growth

External exhaust keeps the mattress surface cooler and the pump quieter

HeelSense™ Technology helps provide further therapy and comfort by decreasing pressure in patient's vulnerable heel area

Request a Quote or Place an Order Today!



**Contact your Inside
Sales Representative**

or



**Contact Customer Service
at 800.333.6900**

Invacare Corporation
www.invacare.com







USA

One Invacare Way
Elyria, Ohio
44035-4190
(800) 333-6900

Canada

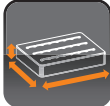570 Matheson Blvd. E.,
Unit 8
Mississauga, Ontario
L4Z 4G4 Canada
(800) 668-5324

microAIR MA800 Alternating Pressure Low Air Loss Mattress System

Model #	 Dimensions	 Weight Capacity	 Flow Rate (Pump)	 Therapy Type	 Alarms	 Limited Warranty
MA800	80", 84" L 36", 42", 48" W 10" H	36" 600 lbs. 42", 48" 1000 lbs.	1275 LPM	Alternating Pressure Pulsation True Low Air Loss	Audio & Visual	1 Year (Cover/Mattress) 1 Year (Power Unit)
MA800	microAIR MA800 36" Alternating Pressure Low Air Loss Mattress System					
MA800B42	microAIR MA800 42" Alternating Pressure Low Air Loss Mattress System					
MA800B4284	microAIR MA800 42"W x 84"L Alternating Pressure Low Air Loss Mattress System					
MA800B48	microAIR MA800 48" Alternating Pressure Low Air Loss Mattress System					
MA800M	microAIR MA800 36" Mattress only					
MA800M42	microAIR MA800 42" Mattress only					
MA800M4284	microAIR MA800 42"W x 84" L Mattress only					
MA800M48	microAIR MA800 48" Mattress only					
MA800P	microAIR MA800 Pump only					

©2019

Invacare Corporation.
All rights reserved.
Trademarks are identified
by the symbols ™ and
®. All trademarks are
owned by or licensed
to Invacare Corporation
unless otherwise noted.
Specifications are subject
to change without
notification.
Form No. 19-026
190205

