

Connect with a NewLeaf Home Medical Advisor at 888.739.4330 for personal assistance  
or to order your very own Invacare micorAIR mattress today!



Shop at [NewLeafHomeMedical.com](https://www.NewLeafHomeMedical.com)

INVACARE® MICROAIR® MA900  
THERAPEUTIC SUPPORT SURFACE



**microAIR®**

# Invacare® microAIR® MA900

Lateral Rotation Low Air Loss Mattress System

**Medium - High Risk  
Stage I-IV**

Alternating Lateral Rotation and True Low Air Loss help provide Pressure Relief and Moisture Management, with built-in side bolsters for added safety



**Highly vapor permeable and oversized pliable quilted nylon top cover** provides low shear, friction and moisture protection



**Cable management system** keeps cords away from the ground to help prevent trip hazards



**Nylon base cover** help to reduce shear and friction during bed articulation







# Features & Benefits

**10 adjustable comfort settings** depending on the patient's therapy requirement

**Auto Firm Mode** provides quick maximum inflation within seconds to aid with transfers and nursing procedures

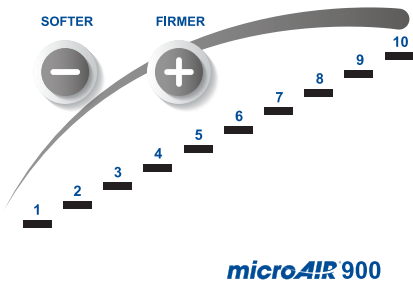

Features **continuous 20° and 40° lateral rotation therapy**, which can help lower the risk of infection, pneumonia and other pulmonary complications in at-risk patients

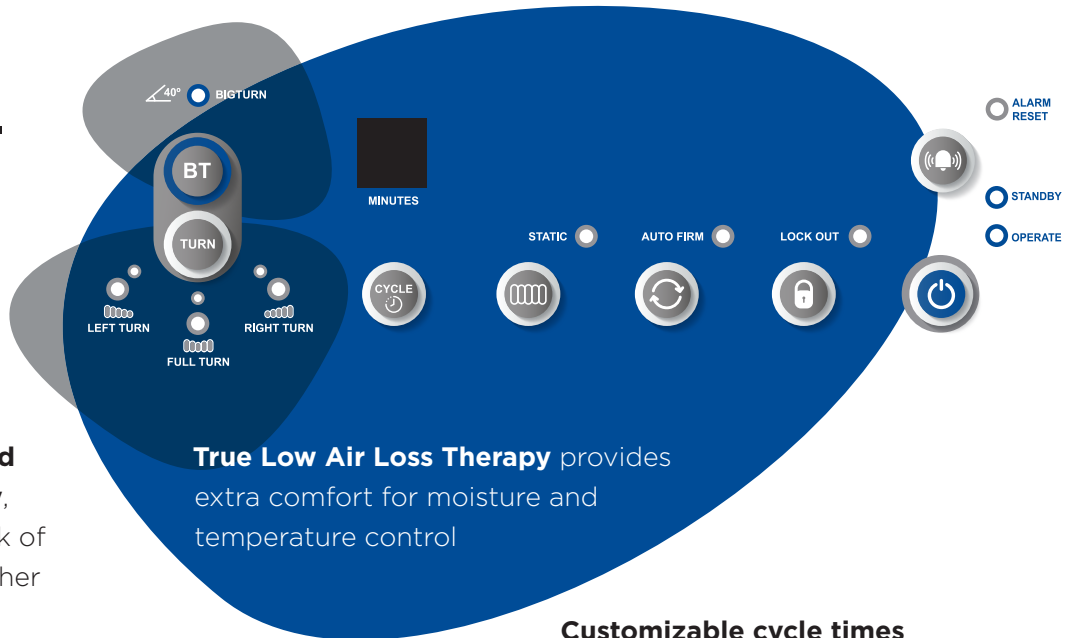

**True Low Air Loss Therapy** provides extra comfort for moisture and temperature control

**Customizable cycle times** to meet patient's varying needs



**Foam base** helps prevent bottoming out in the event of a power failure

**Soft nylon urethane backed cells** provide easy cleaning and helps minimize the proliferation of bacterial growth

**Inflatable built-in side bolsters** with manual connectors can be removed as needed to aid with ease of ingress and egress

**External exhaust** keeps the mattress surface cooler and the pump quieter

# Request a Quote or Place an Order Today!



**Contact your Inside  
Sales Representative**

or



**Contact Customer Service  
at 800.333.6900**

**Invacare Corporation**  
www.invacare.com







## USA

One Invacare Way  
Elyria, Ohio  
44035-4190  
(800) 333-6900

## Canada

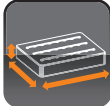570 Matheson Blvd. E.,  
Unit 8  
Mississauga, Ontario  
L4Z 4G4 Canada  
(800) 668-5324

## Invacare microAIR MA900 Lateral Rotation Low Air Loss Mattress System

Model #	 Dimensions	 Weight Capacity	 Flow Rate (Pump)	 Therapy Type	 Alarms	 Limited Warranty
MA900	80" L 36", 42", 48" W 10" H	36" 600 lbs. 42", 48" 1000 lbs.	1275 LPM	Lateral Rotation True Low Air Loss	Audio & Visual	1 Year (Cover/Mattress) 1 Year (Power Unit)
MA900B42	microAIR MA900 42" Lateral Rotation Low Air Loss Mattress System					
MA900B48	microAIR MA900 48" Lateral Rotation Low Air Loss Mattress System					
MA900M	microAIR MA900 36" Mattress only					
MA900M42	microAIR MA900 42" Mattress only					
MA900M48	microAIR MA900 48" Mattress only					
MA900P	microAIR MA900 Pump only					

©2019

Invacare Corporation.  
All rights reserved.  
Trademarks are identified by the symbols ™ and ®. All trademarks are owned by or licensed to Invacare Corporation unless otherwise noted. Specifications are subject to change without notification.  
Form No. 19-028  
190207

