


RHYTHM MULTI

Rhythm Multi Mattress System is a medium to high risk system for patients who are at risk of developing Pressure Ulcers according to your sound clinical judgment. Rhythm Multi can also be used for patients who have an existing stage 1, 2, 3 and 4 Pressure Ulcer, in conjunction with your hospital's policy on pressure area management.

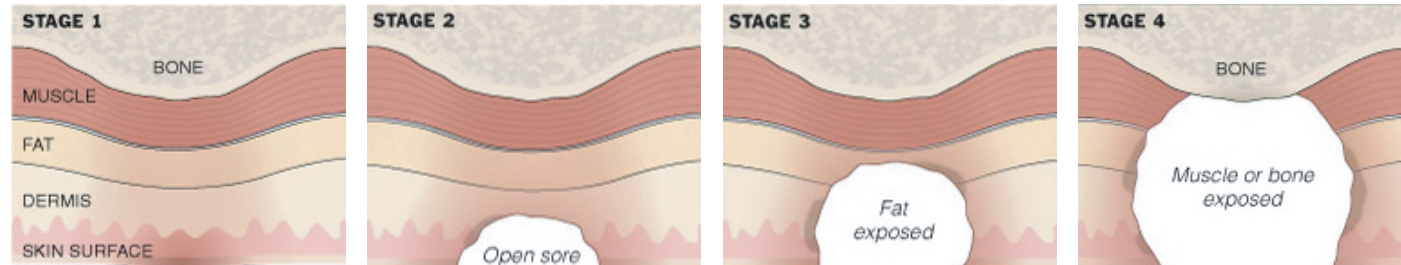

Optional Raised Bolster Cover available for added patient safety. Bolster cover has center access providing easy ingress and egress.

MASTER CONTROL UNIT

Dimension:	(L)17.7"x (W)6.8"x (H)10.8"
Weight:	13.2 lbs
Rating Voltage:	110-120 volts
Max Current:	5A
Rating Frequency:	60Hz
Cycle Time:	3~95 mins
	3 in 1 alternation
Min/Max Pressure:	10~61 mmHg +/- 3~4 mmHg
Protection Type:	Class I, Type BF
	Power outage visual and audio alarm
	Static Function / Pulsation
	Auto-firm for fast mattress inflation
	Keypad lock-out function
Air flow rate:	> 1275 liters per minute

MATTRESS REPLACEMENT

Model No:	Rhythm Multi
Dimensions:	(L)80"x (W)36"x (H)10"
Bariatric Size:	42, 48 inch width
Weight:	28.6 lbs
Weight Capacity:	600, 800, 1000 lbs
Cell Material:	Black Nylon with PU coating
Cell Number:	18 cells
Cover Material:	Blue Nylon with PU coating finish
Base Material:	Black Polyester with PVC coating
	CPR Quick Release System
	Cable Management System on bottom cover
HCPSC Code	E0277

RHYTHM MULTI

Provides you, the caregiver, with peace of mind, and your patient with world class care and comfort




HAPPY HEEL™ Heel Relief
Patent Pending

Perfect for patients and caregivers




HAPPYHEEL™ provides extra comfort settings for the heels

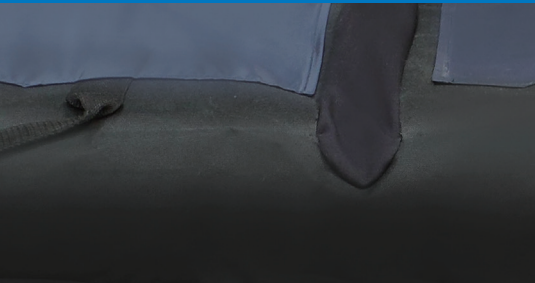

Low slip nylon under mattress base cover to reduce shearing and friction between mattress and bed frame



Nylon woven fabric with PU coating finish

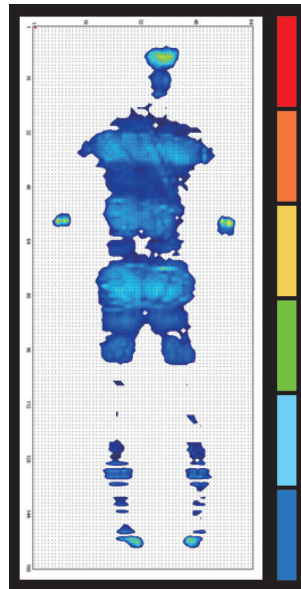


PRESSURE MAPPING

Rhythm Multi provides constant low pressure by redistributing the patient's pressure over a greater area than a standard mattress.

RHYTHM MULTI

- The Rhythm Multi replacement mattress can accommodate patient weights up to 1000 lbs
- PU cell design provides easy cleaning and minimizing the proliferation of bacteria
- 30 seconds of pulsation simulates a massage to assist in maximizing a patient's comfort
- Modular structure design to ensures ease of use and cleaning
- True low air loss provides extra comfort/for moisture and temperature control

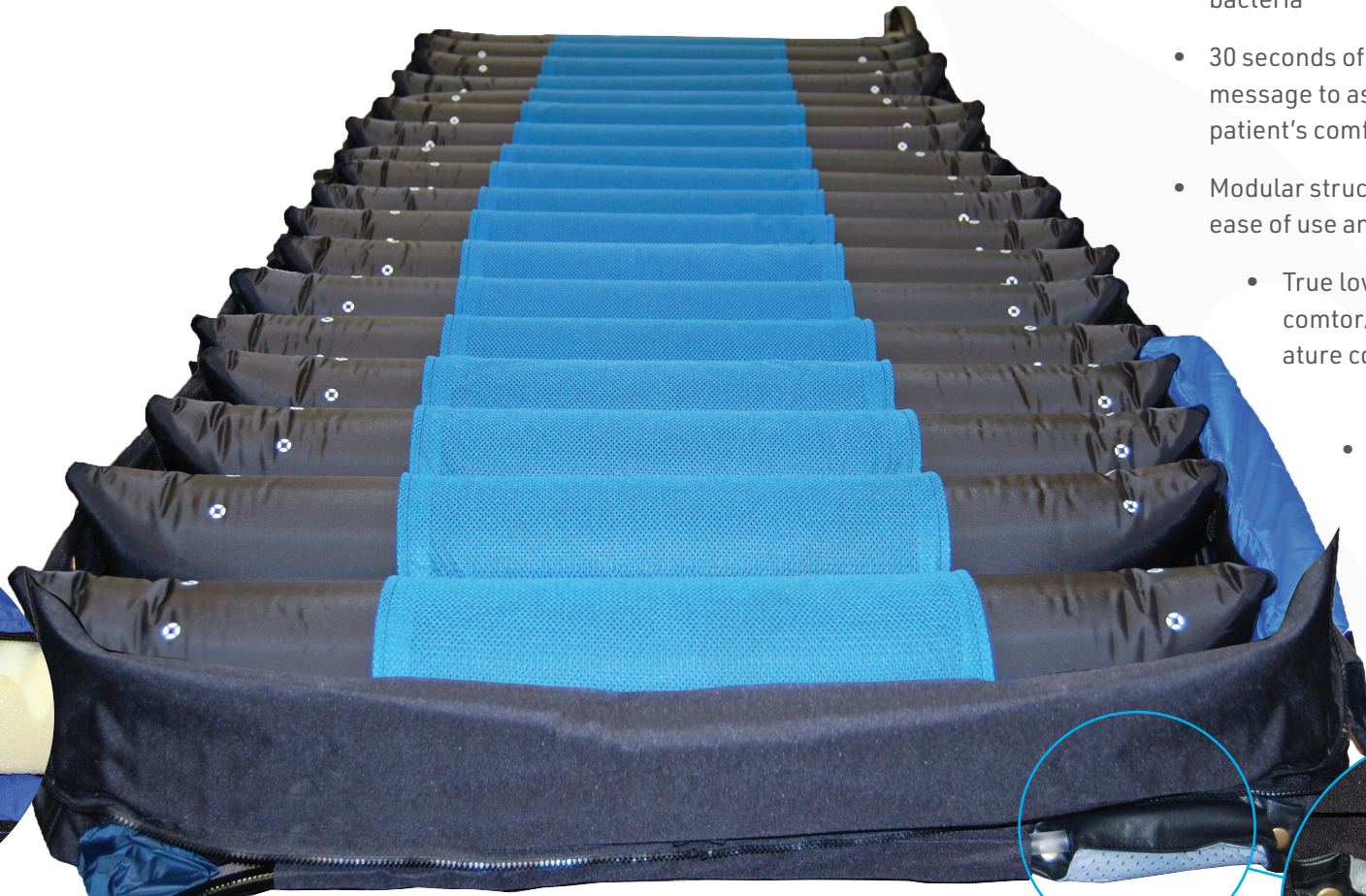

Foam Safety Mat to prevent bottoming out in case of power loss




Exhaust port removes excess heat from pump lowering mattress surface temperature



- provides 3 selectable softness settings in the last 5 air cells to provide protection for patient's fragile heels




Cable management keeps electrical cables off the floor

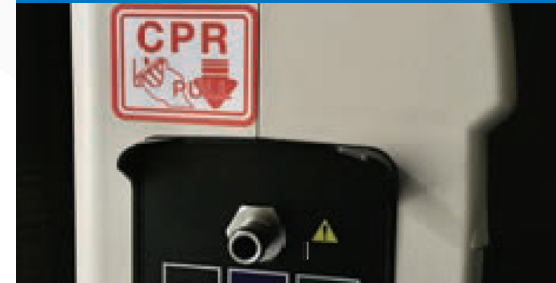

CPR quick release system deflates the system within 20 seconds

10 adjustable comfort settings provides different levels of comfort depending on the patient's therapy required

3 different modes of clinical therapy static/alternation/pulsation

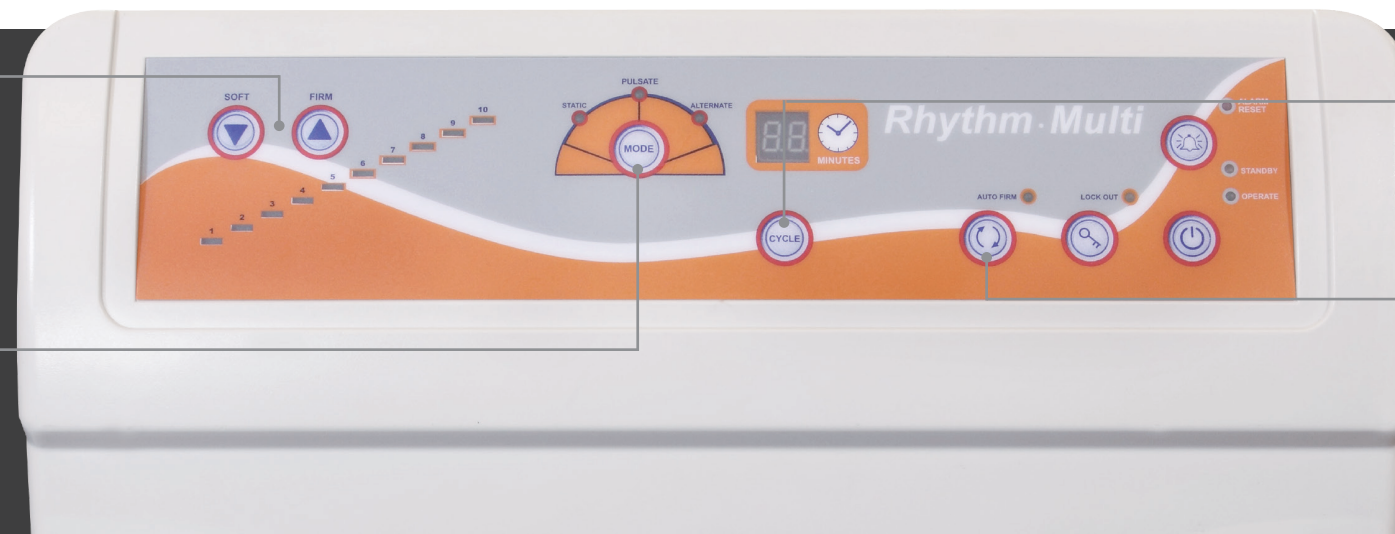

Alternation cycle times from 3-95 minutes

Auto Firm/Patient Care control



Hanging bracket fits on all type of beds

